

# Integra®

Camino® Intracranial Pressure Monitor

Because we are committed to limiting uncertainty, the new Camino® ICP Monitor provides you with an advanced multimodality neuro-monitoring system.



INTEGRA®  
LIMIT UNCERTAINTY

# New CAMINO® monitor

## Advanced ICP System

The Camino® ICP Monitor is the new generation in advanced monitoring featuring an enhanced user interface, ease-of-use, and improved handling.

### Parameters displayed by Camino® monitor

- Intracranial Pressure (ICP): Mean ICP values and continuous waveform
- Intracranial Temperature (ICT): Values

### Compatible with two ICP catheter's technology

- Camino® Fiberoptic Technology
- Camino® Flex Strain Gauge

### Large Touch-Screen Interface

- Intuitive, easy to use menu system
- Widescreen display (7.0 inches):  
Critical numbers and waveform are highly visible

### New Improved Look

- Compact size, monitor weight 3kg
- Ergonomic carrying handle

### Monitor Alarms & Trend Data

- 5 days trend data history of mean ICP values
- Trend data can be displayed at 3, 12, 24, 48 & 120 hrs
- High ICP alarm limit: by numerical display & coloured line on the screen

### Simplified Synchronization

- Single-touch operation
- Simple cal step procedure to ICU bedside monitor



### USB Output (back)

- Convenient and rapid data export through a USB output makes reviewing five days of patient data simple.




Camino® Fiberoptic Technology




Camino® Flex Strain Gauge Technology



### Main Menu

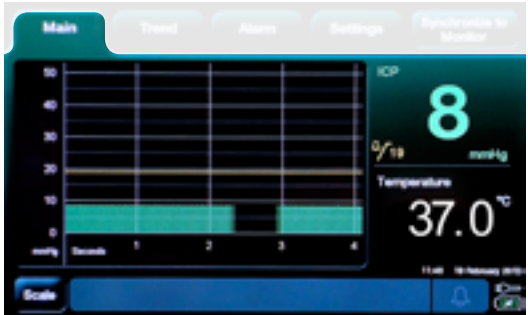


### Trend

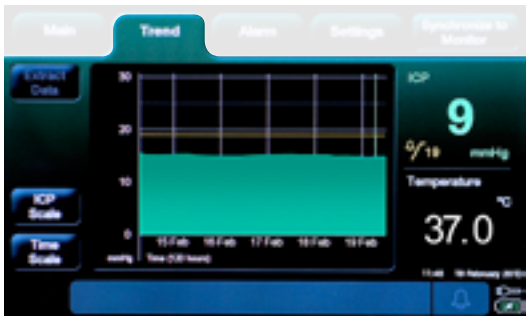


### Alarm



### Settings

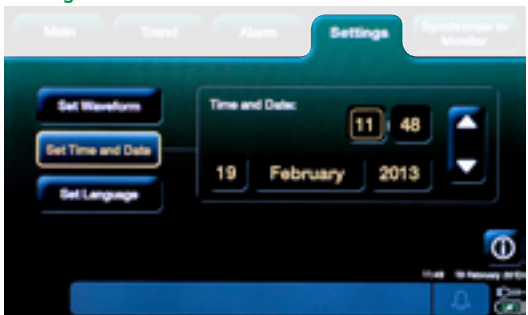

### Zeroing

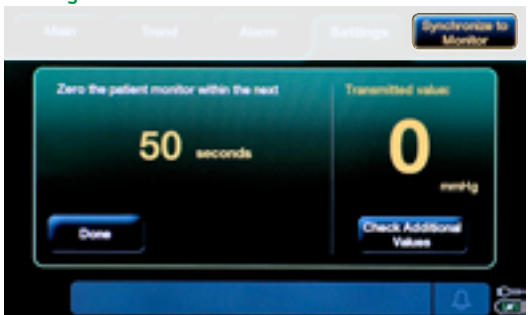

# Integra®

## Camino® Intracranial Pressure Monitor

### Description

The Camino® Advanced Monitor is a compact, portable device for use with Camino® 110-4 series of Pressure/Temperature and Pressure Transducer-Tipped Catheters.

### System includes

- **CAMCABL** Cable (a)  
Extension cable that connects the Fiber Optic Catheters (110-4 series) to the Integra® Camino® ICP Monitor
- **FLEXEXT** Cable (b)  
Extension cable that connects the Integra® Camino® Flex Catheter to the Integra® Camino® ICP Monitor
- **PMIOMPM1** Cable  
Main cable for connecting Integra® Camino® ICP monitor to patient bedside monitor
- **BAT1001** Battery
- **EXPORTCAB** Cable (c)  
USB-to-RS232 adapter cable
- **MONPWR** Power cable



(a)



(b)



(c)

### Product Reference

References	Description
CAM02	ICP Monitor

Compact, portable neuromonitor device that is used for monitoring intracranial pressure (ICP) and intracranial temperature (ICT).

### Indications For Use

- The Integra® Camino® ICP Monitor is indicated for use by qualified neurosurgeons or neurointensivists for measurement of intracranial pressure and temperature.

### Contraindications

- The Integra® Camino® ICP Monitor and its accessories are contraindicated for use in a Magnetic Resonance (MR) environment.

Integra LifeSciences Services (France) SAS  
Sales & Marketing EMEA  
Immeuble Séquoia 2 • 97 allée Alexandre Borodine  
Parc technologique de la Porte des Alpes  
69800 Saint Priest • FRANCE  
Phone: +33 (0)4 37 47 59 00 • Fax: +33 (0)4 37 47 59 99  
emea.info@integralife.com • [integralife.eu](http://integralife.eu)

Customer Services  
International: +33 (0)4 37 47 59 50 • +33 (0)4 37 47 59 25 (Fax) • [cs@integralife.com](mailto:cs@integralife.com)  
France: +33 (0)4 37 47 59 10 • +33 (0)4 37 47 59 29 (Fax) • [custsvcf@integralife.com](mailto:custsvcf@integralife.com)  
United Kingdom: +44 (0)1 264 345 780 • +44 (0)1 264 363 782 (Fax) • [custsvcs.uk@integralife.com](mailto:custsvcs.uk@integralife.com)  
Germany: +49 (0) 2102 5535 6200 • +49 (0)2 102 5536 636 (Fax) • [custsvcg@integralife.com](mailto:custsvcg@integralife.com)  
Benelux: +32 (0)2 257 4130 • +32 (0)2 253 2466 (Fax) • [custsvcbenelux@integralife.com](mailto:custsvcbenelux@integralife.com)  
Switzerland: +41 (0)2 27 21 23 30 • +41 (0)2 27 21 23 99 (Fax) • [custsvcsuisse@integralife.com](mailto:custsvcsuisse@integralife.com)

Integra LifeSciences (Ireland) Limited  
IDA Business and Technology Park  
Sragh, Tullamore, County Offaly • Ireland

©2013 Integra LifeSciences Corporation. All rights reserved. 0032295  
PRODUCTS FOR SALE IN EUROPE, MIDDLE-EAST AND AFRICA ONLY

Distributed by

**INTEGRA**  
LIMIT UNCERTAINTY

Availability of these products might vary from a given country or region to another, as a result of specific local regulatory approval or clearance requirements for sale in such country or region. • Read carefully the Instructions for Use of the product. • Non contractual document. The manufacturer reserves the right, without prior notice, to modify the products in order to improve their quality. • WARNING: Applicable laws restrict these products to sale by or on the order of a physician. • Camino, Integra and the Integra logo are registered trademarks of Integra LifeSciences Corporation in the United States and/or other countries. • All the references numbers mentioned on this document are CE marked according to European council directive 93/42/EEC on medical devices and its relatives, unless specifically identified as "NOT CE MARKED". Class IIb.

CE  
0086  
BSI